MRI Safe Electrodes for Spinal Cord Stimulation (SCS)

Dirk Hartung ^a, Jonathan Scott ^a, Steven McCabe ^a, Peter Single ^b

Introduction

A lead attached to an implant can act as a resonator in an MRI machine. In previous work, we have shown that the worst-case heating of a single-core lead results from strong E-field at the distal end for certain particular lengths. Figure 1 compares simulated worst-case exposure with measured data points obtained in a phantom torso placed in a 3-Tesla MRI machine as function of lead length. We present here further evidence [9] that multi-filar leads spread energy captured from RF fields between filars. We describe how decoy filars and decoy electrodes connected to decoy filars can convey the majority of intercepted energy into human tissue in a distributed fashion.

Methods

Temperature rise around a straight, single filar (Figure 1) has been computed using COMSOL Multiphysics 4.4 and simulation results have been validated by comparison with measurements taken in a 3-Tesla MRI machine (Figure 1). Worst case local SAR values for a twelve filar (electrode) spinal cord stimulator lead in a saline solution (Base Case) have been computed using ANSYS EM Design Suite (HFSS). Then, decoy filars and decoy electrodes have been added to the Base Case (Base Case plus Decoy). The simulation results are presented in Figures 2 and 3.

Results

A single filar lead exposed to an RF field at 128 MHz has an exposure-vs-length function that peaks at a length of about 200 mm. Multi-filar leads show worst-case local SAR values around one particular electrode (Figure 2). Adding decoy filars and non-therapeutic decoy electrodes to therapeutic filars can divert and dissipate a significant fraction of captured RF energy in a distributed way along the lead (Figure 3). Simulated worst case local SAR values have been reduced by a factor of ten for lead lengths ranging from 200 mm to 1200 mm.

Discussion

RF energy is coupled between filars within multi-filar leads. Filars show critical lengths where they capture the majority of the RF energy and may present a risk. Adding decoy filars to therapeutic filars and non-therapeutic decoy electrodes may deflect enough energy to deliver safe operation.

Conclusions

These results promise the possibility of designing multi-wire leads that are unconditionally safe in whole-body MRI scanners. This may require appropriate decoy filars, non-therapeutic decoy electrodes and a distributed release of captured RF energy into human tissue.

Acknowledgement

The authors would like to thank Saluda Medical for funding Waikato University to provide a studentship for Dirk Hartung, and for support in kind to the School of Engineering. The technology is licensed to Saluda Medical.

Figures

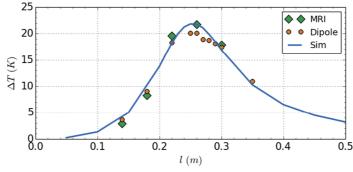


Figure 1 - Simulated and measured temperature rise T near the distal electrode of an insulated wire after 5 minutes in a 3T MRI machine. Measurements from the dipole test method yield a similar profile.

^a University of Waikato, New Zealand

^b Saluda Medical, Australia

Maximum local SAR in saline along lead

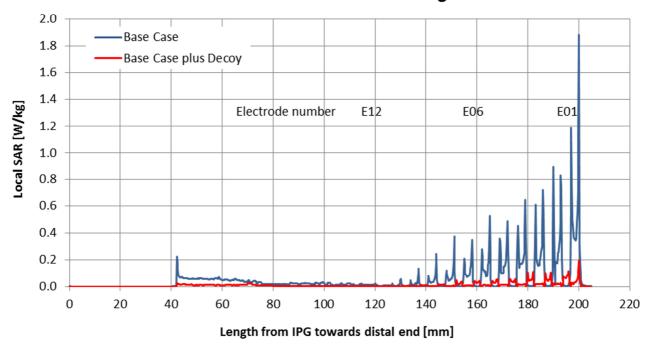


Figure 2 - Local SAR in saline solution normalised to 1 V/m of the exciting electrical field for a spinal cord stimulator with 12 electrodes at the distal end and for a worst case lead length of 200 mm

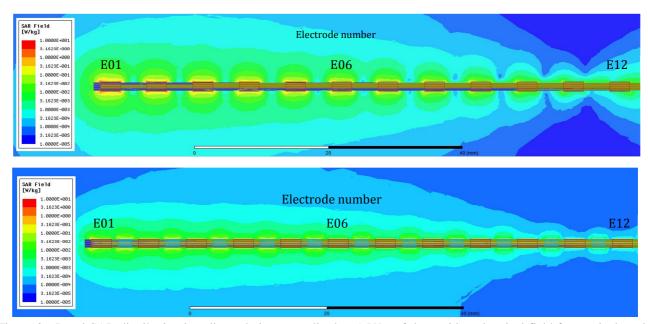


Figure 3 - Local SAR distribution in saline solution normalised to 1 V/m of the exciting electrical field for a spinal cord stimulator with 12 electrodes at the distal end and for a worst case lead length of 200 mm.

Top Base Case

Bottom Base Case plus Decoy

References

- [1] J. A. Nyenhuis, P. Sung-Min, R. Kamondetdacha, A. Amjad, F. G. Shellock and A. R. Rezai "MRI and Implanted Medical Devices: Basic Interactions With an Emphasis on Heating", IEEE Transactions on Device and Materials Reliability, vol. 5, no. 3, pp. 467-480, Sept. 2005.
- [2] R. W. P. King and G. S. Smith, "Antennas in matter: fundamentals, theory, and applications", The MIT Press, Cambridge, MA, 1981.
- [3] F2182-11a Standard Test Method for Measurement of Radio Frequency Induced Heating On or Near Passive Implants During Magnetic Resonance Imaging, ASTM International, West Conshohocken, PA, USA, 2010.
- [4] S. O. McCabe and J. B. Scott, "Cause and amelioration of MRI-induced heating through medical implant lead wires," Proceedings of the 21st Electronics New Zealand Conference, Hamilton, 20-21 Nov. 2014, pp. 34–40.
- [5] S. O. McCabe, "MRI-safe Implantable Electrode Leads", PhD thesis, The University of Waikato, 2016.
- [6] S. O. McCabe and J. B. Scott, "Electromagnetic techniques to minimize the risk of hazardous local heating around medical implant electrodes during MRI scanning," European Microwave Conference 2015, Paris, 7-10 Sept. 2015.
- [7] S. O. McCabe and J. B. Scott, "Implant conductor assembly with improved radio frequency properties", New Zealand provisional patent 714212, Filed on Nov. 16, 2015.
- [8] S. O. McCabe and J. B. Scott, "New MRI-safe implant electrode design," International Microwave Symposium, San Francisco, 22-27 May. 2016, pp. 1-3.
- [9] S. O. McCabe and J. B. Scott, "A Novel Implant Electrode Design Safe in the RF Field of MRI Scanners", IEEE MTT, pp. 1-7, 2017.